



Bariatric Stretchers



Bariatric Stretchers

Designed for large and heavy patients, we've focused on larger, reinforced frameworks, and electronics with increased weight capacity. We also have patient weigh scale and powered-drive options.



Contour Barituff



Contour Orbit-Drive

Contour Barituff

A 500kg SWL extra-wide electric stretcher designed for Bariatric patients.

Features

- 500kg SWL Bariatric stretcher
- 2135 L x 1100mm extra-wide top
- Emergency CPR backrest release
- 15° head and foot down Trendelenburg
- Electric height range: 470 - 850mm
- 980 W x 1900 L x 150mm thick mattress
- 150mm twin wheel castors
- Full body x-ray top option
- Choice of 7 dropside armrest/end rail colours
- MADE IN AUSTRALIA

Product Code


- Barituff - MPT-C-BE
- Barituff-X - MPT-C-BEX

Summary

Having equipment that is robust enough, that you can completely rely on, is the greatest challenge when caring for Bariatric patients. The Contour Barituff has increased componentry strength in the core frame design and operates with higher capacity columns to safely lift 500kg without overload.

The stronger and wider top comfortably caters for large patients, with a choice of mattress styles available to suit the type of departmental use. Stretcher stability is assured through an increased base size. Powered drive and weigh scales can be included to genuinely reduce manual handling risks to the carer.

This stretcher has been built for Bariatric patients.

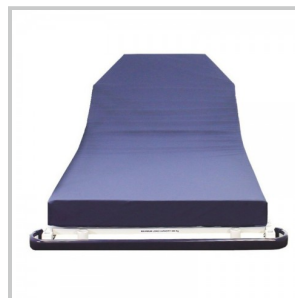

Contour Orbit-Drive

Taking patient transport to a new level in carer welfare with NEW zero effort, 4-way power drive.

Features

- Forward, back and sideways powered drive
- Intuitive controls with emergency stop
- Powerful electric drive centre steer wheel


Product Code

- MPT-CA-EDU

Summary

Take all the human effort out of pushing and jostling stretchers with an intuitive and intelligent drive system designed for easy use by any operator. Natural movement controls the speed and direction of the stretcher with the ability to stop very quickly by pulling back.

Help prevent the most common carer strain injuries by taking the back twist out of moving a stretcher forwards, backwards and sideways. Importantly, we haven't overlooked safety features such as emergency stop and system lock out, which are standard. Add this system to our traditional vertical column stretcher designs such as Contour Barituff, Classic, Deluxe and Multi-X to improve both standard and Bariatric patient transport.

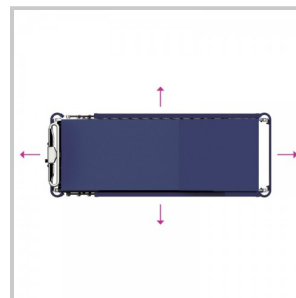

Notes
