



Emergency Stretcher Beds



Emergency Stretcher Beds

In a trauma department, we understand that reliable equipment is critical because it increases patient throughput, improves operating costs and helps keep you in control. We provide emergency stretchers and recovery trolleys that work with your current staff systems and easily adapt to technology upgrades, so you're always covered.

Hospital emergency trolleys

Modsel's exclusive range of hospital emergency trolley designs continues to be a solution for healthcare and nursing professionals.



Contour Portare



Contour Classic



Contour Portare-X



Contour Multi-X



Contour Heli-Transfer

Contour Portare

An electric all-purpose transport and recovery stretcher where Modsel design meets budget constraints.

Features

- Economical transfer trolley
- E-Point steering system
- Emergency CPR backrest release
- 15° head and foot down Trendelenburg
- 300kg safe working load
- 400 - 830mm height range
- Fully electric operation
- 150mm twin wheel castors
- 2135 x 815mm overall size
- 965mm wide top option
- Dropside handrail PVC upgrade option
- Choice of 7 dropside armrest/end rail colours for easy identification
- MADE IN AUSTRALIA

Product Code

- Contour Portare; MPT-C-PE2
- 150mm Wider Top; MPT-CA-WCT
- PVC Dropside upgrade; MPT-CA-CP-PPVC

Summary

Helping patients on and off a general purpose stretcher has never been this easy or safe with a design that allows the mattress base down to 400mm, lower than a standard office chair height!

Transport must always be considered and with our new E-Point system, a strong grip with the floor and solid, perfectly centred steer wheel ensures the stretcher always accurately tracks and turns. We are now able to bring Modsel stretcher design into department areas such as emergency, theatre and day surgery at our most affordable price ever, helping to meet tighter budget constraints and maximising the number of working units available.



Contour Classic

The most flexible custom design platform available for a transport and recovery stretcher.

Features

- Fully customisable stretcher
- A variety of dropside options
- Emergency CPR backrest release
- 15° head and foot down Trendelenburg
- Electric and hydraulic versions available
- 300kg safe working load
- Electric height ranges: from 440 - 810 up to 570 - 940mm
- Hydraulic height ranges: 530 - 760, 565 - 830 or 630 - 960mm
- 200mm castors and Ezi-Trak 5th wheel steering
- 2135 x 815mm overall size
- 705 narrow and 965mm wide top options
- Choice of 7 dropside armrest/end rail colours for easy identification
- MADE IN AUSTRALIA

Product Code

- Electric - MPT-C-CE
- Hydraulic - MPT-C-CH

Summary

The most flexible platform of all our stretcher range. The twin column design allows for several height range choices - we can build in a low height range for patient egress care or higher height ranges for better equipment or surgeon access with less strain during procedures.

Tailor make the Contour Classic with shorter, longer, narrower and wider top sizes and select from a range of dropside styles to suit specialised procedures or the patient type you are working with. To simply put it, there is a high level of customisation possible with many accessory mounting options and solutions.



Contour Portare-X

An electric ultra-low electric trauma stretcher with totally unobstructed x-ray clearance.

Features

- Economical trauma stretcher
- E-Point steering system
- Emergency CPR backrest release
- 15° head and foot down Trendelenburg
- Unique InfineX unobstructed x-ray system
- 300kg safe working load
- 430 - 880mm height range
- Fully electric operation
- 150mm twin wheel castors
- 2135 x 815mm overall size
- 965mm wide top option
- Dropside handrail PVC upgrade option
- Choice of 7 dropside armrest/end rail colours for easy identification
- MADE IN AUSTRALIA


Product Code

- MPT-C-PEX
- 150mm Wider Top; MPT-CA-WCT
- PVC Dropside upgrade; MPT-CA-CP-PPVC

Summary

Helping patients on and off a trauma stretcher is safer and easier when the bed can be lowered down to 430mm, around the standard office chair height! This stretcher features our unique InfineX system which allows totally unobstructed head to toe clearance for effortless and efficient x-ray cassette positioning. Smoothly moving patients helps improve their experience and our E-Point steering delivers a strong grip with the floor with a solid, perfectly centred wheel so the stretcher always accurately tracks and turns.

Emergency departments feature enough trauma. The Modsel x-ray stretchers are designed to be a reliable, cost effective performer, no matter what the shift brings.

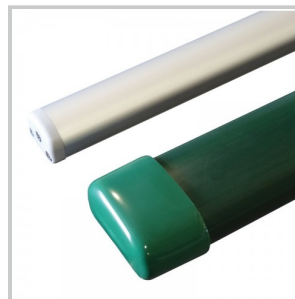

Contour Multi-X

A trauma x-ray stretcher with ability to custom design and build to adapt to any emergency system.

Features

- Fully customisable x-ray trolley
- A variety of dropside options
- Emergency CPR backrest release
- 15° head and foot down Trendelenburg
- Electric and hydraulic versions available
- 300kg safe working load
- Electric height ranges: 545 - 915 or 625 - 995mm
- Hydraulic height ranges: 630 - 860, 665 - 930 or 730 - 1060mm
- 200mm castors and Ezi-Trak 5th wheel steering
- 2135 x 815mm overall size
- 705 narrow and 965mm wide top options
- Choice of 7 dropside armrest/end rail colours for easy identification
- MADE IN AUSTRALIA

Product Code

- Electric - MPT-C-XE-EB
- Hydraulic - MPT-C-XH

Summary

This stretcher features our unique InfneX tracking which allows totally unobstructed head to toe clearance for effortless and efficient X-Ray cassette positioning.

We can help you get the most out of your expensive x-ray systems by building to suit the brand and equipment in use. Tailor make the Contour Multi-X with height range choices and shorter, longer, narrower or wider top sizes. Put safety first and select from a range of dropside styles to suit the department requirements. There is a high level of customisation possible on these Modsel trauma stretchers with many accessory mounting options and solutions



Contour Heli-Transfer

A stretcher to reduce helicopter patient transfer effort outdoors.

Features

- Flat stretcher top with tie down straps and gap bridge
- Fully customised top surface option to suit special purpose
- Can be supplied with or without dropsides
- Standard without backrest, Emergency CPR backrest available
- 15° head and foot down Trendelenburg
- Electric and hydraulic versions available
- 300kg safe working load
- Electric height ranges: 440 - 810 or 530 - 900mm
- Hydraulic height ranges: 530 - 760, 565 - 830 or 630 - 960mm
- 200mm castors and Ezi-Trak 5th wheel steering
- 2135 x 815mm overall size with 705mm narrow and 965mm wide top options
- Choice of 7 dropside armrest/end rail colours for easy identification
- MADE IN AUSTRALIA

Product Code

- Electric - MPT-C-HE
- Hydraulic - MPT-C-HH

Summary

Safely and easily move a helicopter stretcher across onto this transfer stretcher, a platform designed to minimise effort whilst transporting patients over rougher outdoor surfaces. After benefiting from the smooth, flat bridge that lines up with the helicopter door, rapidly fasten the helicopter stretcher securely down with the four quick pull straps.

Reduce the trauma to patients and babies in cribs and the workload for carer's whilst transferring via the large diameter castors and very stable stretcher.



Notes

Notes

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....



Notes
