


Premium Stretcher Mattress

High performance mattress with outstanding durability and patient support within demanding acute environments.

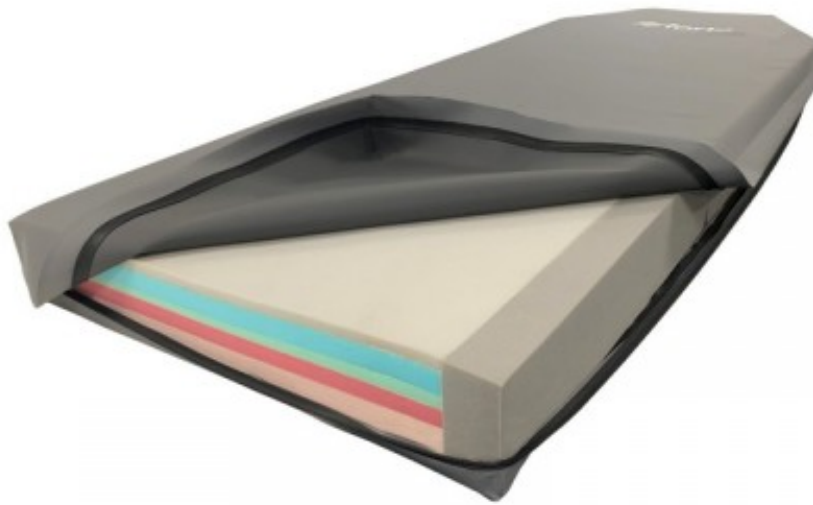

Key Benefits

- **Crucial support:** 5-layered high performance core for demanding environments.
- **Structural stability:** Integrated profile stabilizing sides and unique hinging system.
- **Fluid control:** Double coated, welded polyurethane cover with fluid proof zip.
- **Body distribution:** Superior 110mm mattress thickness with up to 350kg mass range.
- **Peace of mind:** Australian made plus offering a 5 year core foam and 4 year cover warranty.

v-17072023

Premium Stretcher Mattress



Stabilizing sides

High density bars down each side help to prevent roll off and assist the patient during egress by increased side strength of the mattress.




Premium pressure relief

The special combination of foam densities and layers provide superior patient support.



Profile cut sides

Strategically placed reliefs allow this thicker and denser mattress to conform to the stretcher surface contouring.

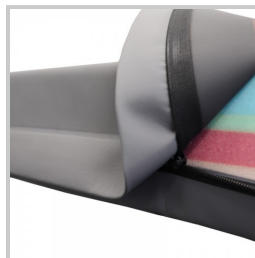

Fluid proof zip

The unique zip design seals when closed to resist fluid penetration, with the added protection of a premium waterfall flap.



Non-slip base

A specially coated base fabric helps to prevent unwanted movement of the mattress and grip the stretcher surface without requiring mechanical fastenings. This allows quick access and fast cleaning of the patient surfaces.

Why Buy

The **Premia Forte** pressure relieving radiolucent support surface is our high spec version of the standard Contour Premia mattress. The specialised core structure provides best protection to patients at risk of injury without compromising support in crucial areas. A Premiflex stretcher cover provides ultimate patient comfort and protection. This mattress delivers an effective outcome - combining both high risk pressure care and longevity for demanding stretcher applications.

Product Code

- Truncated End (hydraulic stretcher); MPT-CM-110-PRF
- Square Ends (electric stretcher); MPT-CM-110-PRPF

Mattress Size

- Standard: 680 W x 1900mm L
- Narrow: 570mm W x 1900mm L
- Wide: 830mm W x 1900mm L

v-17072023