



Bariatric Surgery



Bariatric Surgery

To cater adequately for large and heavy patients, our design team has focused on larger, reinforced frameworks, electronics with increased weight capacity and specialised systems.



Contour Recline Barituff



Oxy Holder



Chair IV Pole



Foot Board



Arm Board



Padded Armrests



Color Options

Contour Recline Barituff

An '1100Lb SWL' Bariatric stretcher chair with an extra-wide top and deeper seat.



Key Benefits

- Multi-purpose treatment and transfer chair for ER, OR surgery, treatment rooms and medical clinics
- 18.9 - 34.6 inch height range with lie-flat and forward / reverse tilt
- 33.5 inch wide x 23.6 inch deep seat designed for large patients
- 38.6 x 50 inch overall (68.9 inch lying flat)
- 4 inch thick premium, removable mattress cushion
- User friendly, long life battery operation
- 1100Lb weight capacity
- Glide-away dropsides with no protrusions
- Steer/brake pedals on both sides
- Sockets on all four corners for many accessories
- Twin wheel 6 inch castors with double steer
- Choice of armrest colours for easy identification

Contour Recline Barituff



Low and safe

Excellent 18.9 inch low height with fold-under drop-sides.

The ultra-low height is combined with drop-sides that glide down vertically under the chair top to provide safe access with no protrusions to catch on the patients legs or clothes.



Easy egress

Forward tilt facilitates getting on and off.

For aged or disabled patients, a gentle incline and the long drop-side armrests may assist in safer transfer.



Lies flat

The backrest and leg rest move independently or together from chair to bed.

You now have the flexibility of sitting a patient up to almost 90 degrees or lying down completely flat.



Head and foot tilt

A generous 15 degrees top tilt both directions.

Position a patient how you need them or quickly lower their head down (and feet up) in an emergency.



Smart transfer

Conveniently stored handles unfold for staff friendly transportation.

Lock in the double steer system, unfold the push handles and move forward effortlessly on the large castors. The unique handle design correctly positions the carer while preventing unwanted forces being applied directly to the patient.



Large top design

The Barituff has a 7.8 inches wider x 5.5 inches deeper seat than the standard Recline.

The seat size is increased to 33.5 inch wide x 23.6 inch deep, easily accommodating larger patients. Use our range of colour coding choices if required to assist in quick identification.



Clear handset

Large simple functions with colour assistance.

Position the patient quickly and easily with the touch of a button - either by independently moving sections or applying one touch contouring.

Why Buy

Having equipment that is robust enough, that you can completely rely on, is the greatest challenge when caring for Bariatric patients. We have increased componentry strength in the core frame design which operates with higher capacity columns to safely lift 1100Lb without overload. A wider top with a deeper seat has been standardized to comfortably cater for large patients. Heavy duty, dark coloured vinyl is used on the mattress cushions to combat wear and staining.

This stretcher chair has been built for Bariatric patients.

Product Code

- MPC-C-REB

Oxy Holder

Stainless steel C size oxygen bottle holder which can be mounted to all four corners of the chair.

Key Benefits

- Powder coated stainless steel construction
- Holds C-size oxygen bottle
- Low, high and backrest mount positions

Product Codes

- Vertical Corner Mount - MPC-CA-VO2
- Vertical Corner Upright Mount - MPC-CA-VUO2
- Vertical Oxy Backrest Mount - MPC-CA-BO2

Why Buy

To cater for varied departmental requirements for specific oxy bottle positioning, we have developed a range of options. Some prefer fixed corner mount positions and others require a swinging style mounted from the backrest. There is also the option to support bottles vertically either above or below the patient surface.



Chair IV Pole

Height adjustable, non-rotational pole that fits all four corners of the chair and with under-top storage option.

Key Benefits

- A strong and safe hook design
- Light and easy 'lift out' auto-lock hanger
- Non-rotational design to prevent tube tangling
- High strength locking for heavy loads
- Tough construction from stainless steel and nylon
- Easy 'pull-down' release from any angle
- Large diameter stand body for ease of pump mounting
- Friction retainer mount for effortless use
- Under top storage option available
- 27.6 inches up to 44.7 inches total height
- 22 lbs weight capacity, 11 lbs per hook
- Only 2.2 lbs total weight

Product Codes

- IV Pole; MPC-CA-IVP
- IV Storage Holder; MPC-CA-IVH

Why Buy

Constructed tough from stainless steel and nylon, these clever non-rotational IV poles are designed to withstand the harshest use and make daily life easier by helping to prevent the tangling of critical soft tubing and cords. Simply lift out the hanger to any height and it automatically locks into an infinite position. Pull the lock body down from any direction to release and lower. This removable IV fits into all four corners of a chair top with the option to store underneath if there is provision.



Versatile



IV storage option

Foot Board

A removable board that locks into a padded foot rest or folds out for extra length as a leg support.

Key Benefits


- Covered in heavy duty vinyl
- Fits standard Recline leg support
- Use as a foot rest or fold down to lengthen the leg support

Product Code

- MPC-CA-DDF

Why Buy

Clip in this foot board to provide a comfortable padded rest for the patients feet to help restrain them whilst moving the chair. It is covered in hard wearing vinyl that is easy to clean. The board can be folded down if the legrest is raised to provide extra length and support for the patient legs.

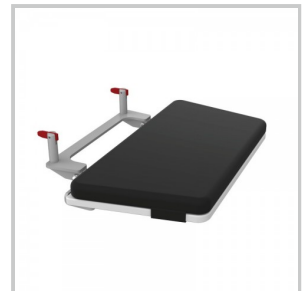

Arm Board

Lightweight armboard that clips onto the dropside hand rail and is padded for patients comfort.

Key Benefits

- 25mm high density foam covered in PU material
- Solid medical grade plastic base board
- Clips on both inside and outside of the dropside
- Adjustable angles for best patient positioning

Product Code

- MPC-CA-SAB

Why Buy

Simply and securely lock the arm board in a choice of locations ranging from the front to the back of the dropside that suits the size and position of the patient. It may also be locked into a variety of angles to help the carer's access. The design is light to carry and store, quick to use and easy to clean.



Padded Armrests

Thick clip-on pads to soften the dropside armrests.

Key Benefits

- High density foam core with slow recovery outer
- Covered in welded, pressure relieving stretch material
- 14.5 long x 4.9 wide x 8.3 inches high
- 1.6 inches of padding on top

Product Code

- MPC-CA-PDA

Why Buy

The Contour Recline Chair has been designed as a practical, clinical product that is functional and easy to clean. In some cases a patient may require extra support or comfort if spending longer periods of time in the chair. These easy to handle and store pads clip on over the dropside to widen and soften the armrest surfaces.



Color Options

A selection of colors to customize the Contour Recline procedure chairs for easy identification.



Key Benefits

- Tough extruded covers for a comfortable hand grip and bump protection.
- Used on chair dropside hand rails.
- Easy to replace for colour change or wear and tear.

Color Options



Flame Red



Midnight Blue



Deep Violet (Purple)

Product Codes

- Flame Red - MPT-CC-FR
- Midnight Blue - MPT-CC-MB
- Deep Violet (Purple) - MPT-CC-DV

Notes

Notes

Notes
