

Contour E-Vertex

A fully electric stretcher with effortless head positioning and comfortable access.



Key Benefits

Full flexibility: The infinitely adjustable headrest up and down tilt and excellent raise height allows the surgeon to place a patient in exactly the right position.

Safety and support: The compact, rigid locking headrest with increased leg clearance brings patients close and keeps them stable while minimizing surgeon neck and back strain.

Comfort and control: Instantly increase patient comfort and enhance access with one-touch, electronic control of the scalloped or flat cushioned headrest; and lightweight removable frame with cushions to support a pillow during recovery.

Slim and clean: The ultra-thin design hides all mechanical operation and has a one-piece, gap-free top surface for ease of cleaning and infection control.

Application

Ideal for Eye, ENT, plastics and dental care.

Why Buy?

We know precision care requires 100% confidence in the equipment being used. That's why the Contour E-Vertex combines state-of-the-art technology and ergonomic design to deliver exceptional support for the performance of procedures requiring close patient proximity, control for the practitioner and patient comfort. The slimline design of this extremely stable, rigid locking headrest provides excellent leg clearance – bringing the patient closer. While the touch button electronic positioning system* allows a surgeon or their support team to quickly, safely and precisely position a patient's head. [Also available as a procedure chair, see the Contour Recline E-Vertex here.](#) *International patents pending

v-19102022

Contour E-Vertex



Slim and clean

The ultra thin design safely hides all mechanical operation and has a one-piece, gap free top surface for easy cleaning.



Effortless control

A clearly marked handset infinitely positions the headrest with the touch of a button.



Raise and lower

Adjust the height up and down while remaining level.




High tilt

Excellent raise height with up and down tilt adjustment.



Low tilt

A practical low height with up and down tilt adjustment.

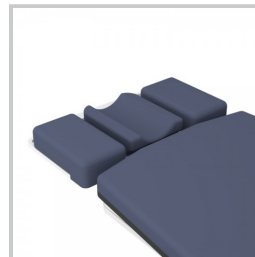

Pillow Support option

Lightweight removable frame with cushions to help support a pillow during recovery.



Wrist Support option

Coated, 11 inches wide, height adjustable and removable rest.



Fold-back design

The wrist support flips away to allow quick and easy access to the patient.

Specifications

Size: 87 x 32 inches overall with 38 inch wide option

Height Range: 20 - 34 inches

SWL: 660Lb

Castor Size: 8 inch castors and Ezi-Track 5th wheel steering

Package Includes: 75mm premia mattress, vertical oxy holder, folding IV, USB Port, 2nd Handset, Emergency CPR Release

Product Codes

Contour E-Vertex; MPT-C-VE-E

Wrist Support; MPT-CA-OWS

Pillow Support; MPT-CA-CO-OSS

Flat Headrest Cushion; MPT-CA-FHC

v-19102022