



Stretchers



Stretchers


Our stretchers and accessories are designed through direct feedback from healthcare professionals and manufactured by us, to ensure you receive a quality product that enables you to deliver the most efficient and effective patient care.



Contour Classic



Stretcher IV Pole

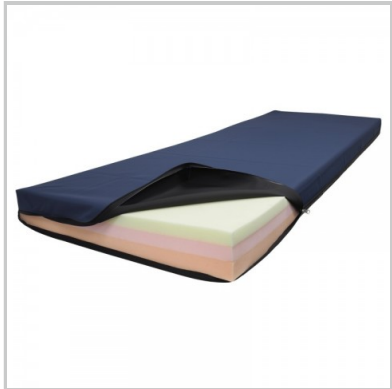

Stretcher Mattress



Oxy Holders



Dropdown Foot Board

Contour Classic

The most flexible custom design platform available for a transport and recovery stretcher.

Features

- Fully customisable stretcher
- A variety of dropside options
- Emergency CPR backrest release
- 15° head and foot down Trendelenburg
- Electric height ranges: from 440 - 810 up to 570 - 940mm
- 660Lb safe working load
- 6 inch castors and Ezi-Trak 5th wheel steering
- 84 inch x 32 inch overall size
- 28 inch narrow and 38 inch wide top options
- Choice of 7 dropside armrest/end rail colors for easy identification


Product Code

- MPT-C-CE

Summary

The twin column design of this stretcher allows for several height range choices - we can build in a low height range for patient egress care or higher height ranges for better equipment or surgeon access with less strain during procedures.

Tailor make the Contour Classic with shorter, longer, narrower and wider top sizes and select from a range of dropside styles to suit specialised procedures or the patient type you are working with. To simply put it, there is a high level of customisation possible with many accessory mounting options and solutions.

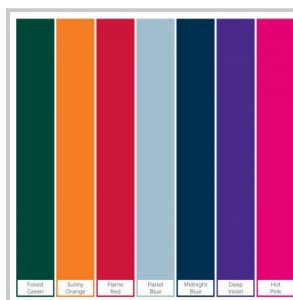

Stretcher IV Pole

Include a standard removable or fold down IV Pole with your stretcher.

Features

- A strong and safe hook design
- Light and easy 'lift out' auto-lock hanger
- Non-rotational design to prevent tube tangling
- High strength locking for heavy loads
- Tough construction from stainless steel and nylon
- Easy 'pull-down' release from any angle
- Large 1.1 inch diameter stand body for ease of pump mounting
- Friction retainer mount for effortless use
- Removable and permanent fold down options
- 28 inch up to 44 inch total height
- 22Lb weight capacity, 11Lb per hook
- Only 2.2Lb total weight

Product Codes

- Standard IV Pole - MPT-CA-IVP

Summary

Constructed tough from stainless steel and nylon, these clever non-rotational IV poles are designed to withstand the harshest use and make daily life easier by helping to prevent the tangling of critical soft tubing and cords. Simply lift out the hanger to any height and it automatically locks into an infinite position. Pull the lock body down from any direction to release and lower. The standard removable IV fits into all four corners of a stretcher top and stores underneath in the base if there is provision. The folding IV fits onto either end of the stretcher top and can be installed to fold down either to the left or right.



Stretcher Mattress

High grade mattress with excellent durability and pressure relieving comfort for patient stretchers.

Features



- A thickness that aid egress and imaging
- 3 effective, layered foams
- Welded 'double coated' stretch cover
- Non-slip base design
- Permanent care label
- New technology fluid proof zip

Product Code

- PT-CM-100-PRP

Summary

The **Contour Premia** mattress has been designed for an increased patient weight capacity and comfort by using three quality foam layers, a high density layer to assist in preventing patients from bottoming and a slow recovery, memory foam upper layer to reduce pressure. It is supplied with a quality welded, stretch and fluid proof cover with permanent label to assist in caring for the mattress.

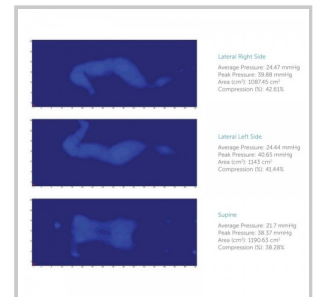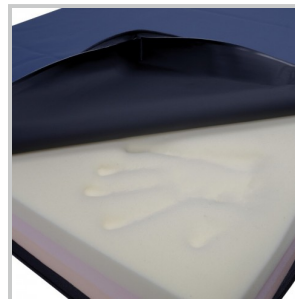

Oxy Holders

A selection of stainless steel oxy bottle holders that can be mounted vertically or horizontally.

Features

- Vertical basket size for up to 4.7 inch diameter bottle; 12.6 deep x 6 inches wide
- Horizontal basket size for up to 6.3 inch diameter bottle; 23.6 long x 7 inches wide
- Horizontal basket can vary in hanging height depending on stretcher type; typically 11 inches

Product Codes

- Vertical Rail Mount - MPT-CA-VO2
- Vertical Corner Mount - MPT-CA-VCO2
- Vertical Corner Upright Mount - MPT-CA-VCUO2
- Horizontal Oxy Holder - MPT-CA-HO2

Summary

Due to so many wide and varied departmental requirements for specific oxy bottle sizes and positioning, we have developed a range of options. Larger bottles can be laid horizontally to ensure they are kept out of the way and smaller bottles are commonly mounted vertically. Some prefer 'fixed', non-moving positions and others require 'swinging' options that allow the bottles to be moved aside if they become an obstruction during transport. There is also the option to support bottles vertically either above or below the patient surface.



Dropdown Foot Board

Padded fold-down foot board to increase the length of the patient surface.

Features

- High density foam covered in welded stretch material
- Tough board construction in a one-piece hospital grade plastic
- Stainless steel lift and lock folding hinge
- 21 wide x 13 inches high (10 inch high cushion)
- Adds an extra 12 inches to the length of a mattress

Product Code

- MPT-CA-DDF

Summary

Padded for patient comfort, the removable dropdown foot board has molded handles for pushing the stretcher and with a simple 'lift and lay' action, folds down to quickly increase the length of the mattress area.



